




BIOHORIZONS®

tapered internal & tapered internal plus

The chart below illustrates the recommended drilling sequence for the BioHorizons Tapered Internal and Tapered Internal Plus Implant systems. Clinical judgement with respect to each individual patient must supersede any recommendation made.

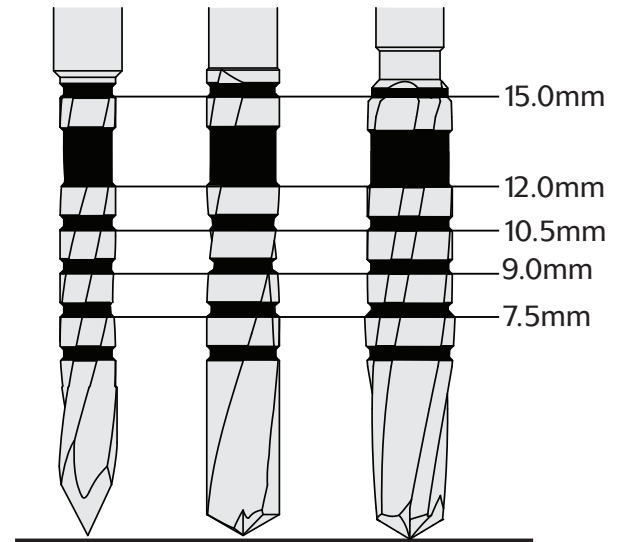

DRILL SEQUENCE

3.0mm implant diameter	3.8mm implant diameter	4.6mm implant diameter	5.8mm implant diameter				
Depth Drills	Width Increasing Drills (implant diameter specific)		Depth Gauge	Crestal Bone Drill (site specific)	Bone Tap (site specific)	Implant Driver	
Use to set osteotomy depth	Incrementally widen the osteotomy without creating excessive heat		Use to verify depth of osteotomy	Use for crestal placement	Tap required in dense bone; stop when resistance gives or desired depth is reached	Place implant matching the length of the prepared osteotomy	

