


# Tapered Implant Site Dilators

When bone quality is at a premium




The Tapered Implant site dilators are used to create or enlarge an osteotomy in soft maxillary bone. These instruments compress the bone laterally rather than removing valuable bone from the surgical site, creating a more dense bone-to-implant interface. The Tapered site dilators match the geometry of the Tapered Surgical Drills and are used sequentially until the desired osteotomy is developed.

- Compatible with Tapered Internal, Plus, 3.0 and Tissue Level
- Fully anodized and laser-marked for easy identification
- Ultra light for easy handling and extended use
- Now includes 2.8mm site dilator

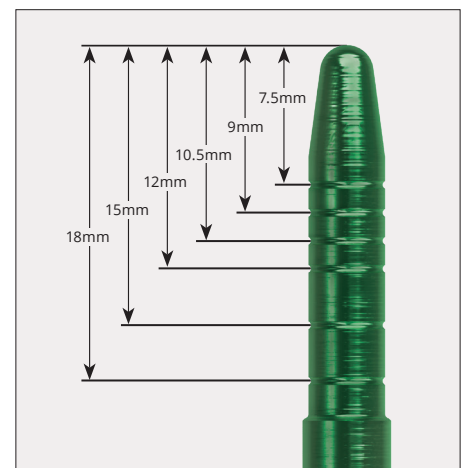

**TODKIT2**  
Tapered Offset  
Site Dilator Kit

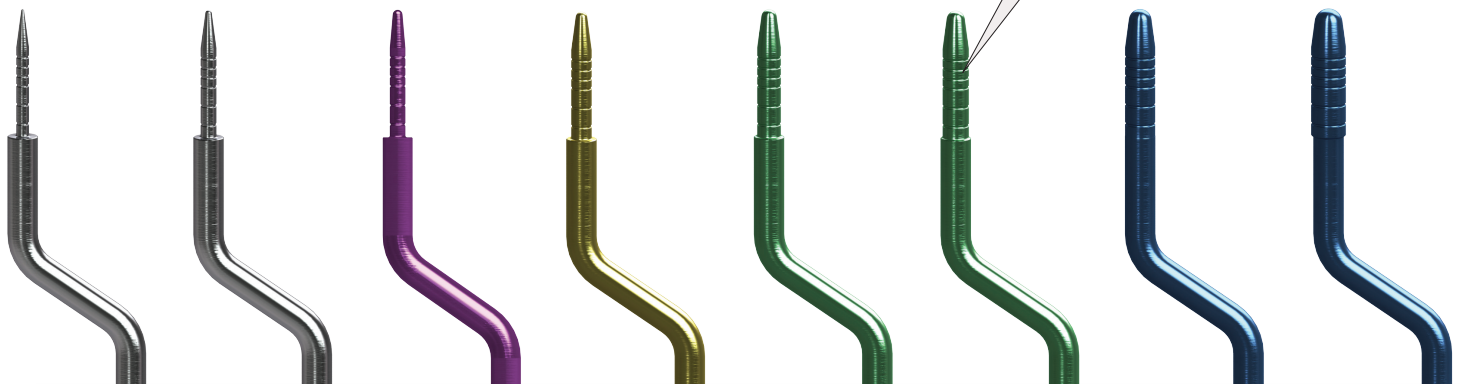

TOD2020	TOD2025	TOD2028	TOD2032	TOD2037	TOD2041	TOD2047	TOD2054
2.0mm Offset Tapered Implant Site Dilator	2.5mm Offset Tapered Implant Site Dilator	2.8mm Offset Tapered Implant Site Dilator	3.2mm Offset Tapered Implant Site Dilator	3.7mm Offset Tapered Implant Site Dilator	4.1mm Offset Tapered Implant Site Dilator	4.7mm Offset Tapered Implant Site Dilator	5.4mm Offset Tapered Implant Site Dilator

BioHorizons® is a registered trademark of BioHorizons. Not all products are available in all countries. ©BioHorizons. All Rights Reserved. ML0227 REV C NOV 2021

customer care: **1.866.468.8338**  
order 24/7 at [castore.biohorizons.com](http://castore.biohorizons.com)

